

LUSCENT-R™ ANCHORS

TRANSLUCENT ENDODONTIC POSTS

Next Generation Fiber Posts with Increased Radiopacity

The Advantages

- Radiopaque
- Superior esthetics
- Dentin-like elasticity
- Biocompatible
- Priced right
- Expedited procedure
- Atraumatic retrieval
- Swedish quality

Without compromise, LUSCENT-R offers all the benefits expected from brand leading, light-transmitting fiber posts, now with higher radiopacity levels for clear & crisp x-rays.

Superior mechanical strength without compromising flexibility

Light transmitting

Textured surface produces strong bond to dentin



STARTER KIT

LUSCENT-R™ ANCHORS

Starter Kit 3, With Standard Conical Reamers

- LUCR-K3** 15 LUSCENT-R Anchors, color coded depth markers
(5 of each size: small, medium, large)
 3 Assorted Standard Conical Reamers (one of each: RUH-4, RUH-5, RUH-6)
 15 Core Forms-to-Fit



REFILLS

LUSCENT-R™ ANCHORS

10 LUSCENT-R Anchors + 10 Assorted Core Forms-to-Fit/pkg

- LUCR-L10** Luscent-R Anchors Large (blue depth markers)
LUCR-M10 Luscent-R Anchors Medium (red depth markers)
LUCR-S10 Luscent-R Anchors Small (yellow depth markers)
LUCR-XS10 Luscent-R Anchors Extra Small (white depth markers)
LUCR-XXS10 Luscent-R Anchors XX-Small (green depth markers)
LUCR-XXXS10 Luscent-R Anchors XXX-Small (purple depth markers)



ACCESSORIES

CORE FORMS-TO-FIT, REFILLS

100 pcs/pkg

- LX-F** Assorted
LX-FS Small
LX-FM Medium
LX-FL Large



LUSCENT-R™ ANCHOR SIZES

XXXS: LENGTH 15 MM & Ø 1,00 MM		S: LENGTH 15 MM & Ø 1,45 MM	
XXS: LENGTH 15 MM & Ø 1,15 MM		M: LENGTH 17 MM & Ø 1,60 MM	
XS: LENGTH 15 MM & Ø 1,30 MM		L: LENGTH 19 MM & Ø 1,75 MM	

DENTATUS LUSCENT ANCHORS®

FIBERGLASS RESIN ANCHORS

Premiere Product for Strong Aesthetic Restorations

The Advantages

- The original fiber anchors from Dentatus
- Excellent esthetics
- One step procedure
- Creates anchor-core support

For strong aesthetic anchorage of dental cohesive restorations. One-Step procedure for simultaneously curing adhesives and composites within the confines of canals creates a mono-bloc strength anchor-core support.

Color coded for fast size indication

Provide extremely high adhesion to composite materials

6 diameters to fit very slim and very large canals



STARTER KIT

LUSCENT ANCHORS

Starter Kit 3, With Standard Conical Reamers

- LUC-K3** 15 Luscent Anchors, color coded depth markers (5 of each size: small, medium, large)
 3 Assorted Standard Conical Reamers (one of each RUH-4, RUH-5, RUH-6)
 1 Probos I Pathfinder #3
 15 Assorted Core Forms-to-Fit



REFILLS

LUSCENT ANCHORS

5 Anchors + 6 Ass. Core Forms-to-Fit/pkg

- LUC-L5** Luscent Anchors Large (blue depth markers)
LUC-M5 Luscent Anchors Medium (red depth markers)
LUC-S5 Luscent Anchors Small (yellow depth markers)
LUC-XS5 Luscent Anchors Extra Small (white depth markers)
LUC-XXS5 Luscent Anchors XX-Small (green depth markers)
LUC-XXXS5 Luscent Anchors XXX-Small (purple depth markers)



15 Anchors + 15 Ass. Core Forms-to-Fit/pkg

- LUC-L15** Luscent Anchors Large (blue depth markers)
LUC-M15 Luscent Anchors Medium (red depth markers)
LUC-S15 Luscent Anchors Small (yellow depth markers)
LUC-XS15 Luscent Anchors Extra Small (white depth markers)
LUC-XXS15 Luscent Anchors XX-Small (green depth markers)
LUC-XXXS15 Luscent Anchors XXX-Small (purple depth markers)



ACCESSORIES

CORE FORMS-TO-FIT, REFILLS

100 pcs/pkg

- LX-F** Assorted
LX-FS Small
LX-FM Medium
LX-FL Large



LUSCENT ANCHOR SIZES

XXXS: LENGTH 15 MM & Ø 1,00 MM



S: LENGTH 15 MM & Ø 1,45 MM



XXS: LENGTH 15 MM & Ø 1,15 MM



M: LENGTH 17 MM & Ø 1,60 MM



XS: LENGTH 15 MM & Ø 1,30 MM



L: LENGTH 19 MM & Ø 1,75 MM



DENTATUS TWIN LUSCENT ANCHORS®

FIBERGLASS RESIN DENTAL ANCHOR WITH AN INNOVATIVE DESIGN

A Winning Combination of Light Transmission, Attractive Esthetics and Twice the Retention

The Advantages

- Design provides extra high retention
- Glass fiber material for high aesthetics
- High adhesion
- Twice the invention
- Twice the retention

The Twin Luscent Anchor is based on the Luscent Anchor. The design is a visible assurance against accidental debonding of adhesive and resin-core materials.

Dentatus has raised the standard of excellence with the Twin Luscent Anchors. This innovative design is your visible assurance against accidental debonding of adhesive and resin-core materials. The slim mid-section creates a "physical choke." The vent groove eliminates air resin entrapment and prevents rotational dislocation. It all adds up to a winning combination of light transmission, attractive esthetics and twice the retention.*

**Al-Harbi & Nathanson, Boston U. Presented at AAED/IFED MTG, 10/3/02*

Color coded for fast size indication

Hour-glass shape provides improved retention

Longitudinal groove provides anti-rotational interlocking



STARTER KIT

TWIN LUSCENT ANCHORS

Starter Kit 1, With Standard Classic Reamers

- LUT-K1** 15 Twin Luscent Anchors, color coded depth markers (5 of each size: small, medium, large)
 3 Assorted Standard Classic Reamers (one of each RUA-4, RUA-5 and RUB-6)
 15 Assorted, Core Forms-to-Fit



REFILLS

TWIN LUSCENT ANCHORS


5 Anchors + 6 Ass. Core Forms-to-Fit/pkg

- LUT-L5** Luscent Anchors Large (blue depth markers)
LUT-M5 Luscent Anchors Medium (red depth markers)
LUT-S5 Luscent Anchors Small (yellow depth markers)
LUT-XS5 Luscent Anchors Extra Small (white depth markers)



15 Anchors + 15 Ass. Core Forms-to-Fit/pkg

- LUT-L15** Luscent Anchors Large (blue depth markers)
LUT-M15 Luscent Anchors Medium (red depth markers)
LUT-S15 Luscent Anchors Small (yellow depth markers)
LUT-XS15 Luscent Anchors Extra Small (white depth markers)

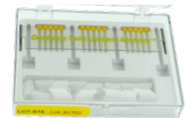

ACCESSORIES

CORE FORMS-TO-FIT, REFILLS

100 pcs/pkg

- LX-F** Assorted
LX-FS Small
LX-FM Medium
LX-FL Large



TWIN LUSCENT ANCHOR SIZES

XS: LENGTH 16,7 MM & Ø 1,26 MM



M: LENGTH 16,7 MM & Ø 1,54 MM



S: LENGTH 16,7 MM & Ø 1,40 MM



L: LENGTH 16,7 MM & Ø 1,68 MM

