

# DENTATUS REAMERS

CREATING THE PERFECT POSTULATE

## The Perfect Match to Dentatus Post Systems

### The Advantages


- For a perfect match between the root and the canal posts
- High quality stainless steel
- Wide choice of sizes

Dentatus reamers are designed to create a perfect match between dental root canals and the shapes and sizes of the corresponding posts. Made in Sweden from surgical high quality steel for high quality performance.

Reliable and long lasting performance

Sizes are indicated by numbers of grooves on the shank

Sharp cutting edges to create perfect seat for Dentatus Posts

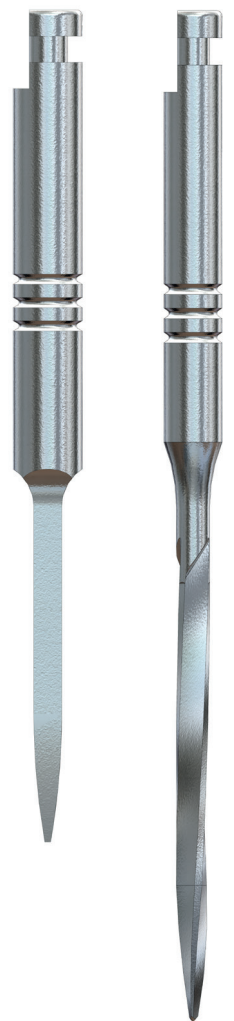

## Helix Reamers

Dentatus Helix Reamers are made of toughened surgical quality stainless steel, the reamers are reliable and long-lasting in their performance, with high cutting capacity, smooth cutting action and improved transport of cutting debris.

They have sharp cutting edges to create perfect seats for Dentatus root canal posts. The Helix reamers are available in both classic and conical shapes and in several diameters and lengths.

Sizes are indicated by number of grooves on the shank. Number of grooves correspond to the same size number for post size on the Dentatus posts.

Use Helix Classic Reamers for Classic Surtex, Twin Luscent and Luminex posts, Helix Conical Reamers for Luscent.



## Standard Reamers

The Original Dentatus reamers. Made in the same high quality stainless steel as all our reamers, for long lasting performance.

Sizes are indicated by number of grooves on the shank, which correspond size number for post size on the Dentatus metal posts.

Use Standard Classic Reamers for Classic Surtex, Twin Luscent and Luminex posts, Standard Conical Reamers for Luscent.



## Special Reamers

**Dentatus Probos I Pathfinder** is used for initial and safe removal of root filling material. Its square cross-cut working tip minimizes the risk for perforation of the dentinal walls. A built-in safety feature makes the Pathfinder to fracture next to the shank, if it gets stuck in the root canal.

**Dentatus Probos II Reamer/Router** has a special profile design which creates additional indents in the dentinal wall, when moved laterally. The indentations add extra surface to the dentinal wall, improving retention. Sizes are indicated by number of grooves on the shank. Number of grooves correspond to the same size number for post size on the Dentatus metal posts.



## Corresponding Posts & Reamers

### Corresponding reamers to your Classic Surtex® posts

HELIX CLASSIC: RSA #1-6, RSB #1-6, RSC #1-5  
STANDARD CLASSIC: RUA #1-6, RUB #1-6.

### Corresponding reamers to your Luscent Anchors®

HELIX CONICAL: RSG #3-6, RSH #3-6, RSI #3-5  
STANDARD CONICAL: RUG #1-6, RUH #4-6.

### Corresponding reamers to your Twin Luscent Anchors®

HELIX CLASSIC: RSA #3-6, RSB #3-6, RSC #3-5  
STANDARD CLASSIC: RUA #3-6, RUB #3-6.

### Corresponding reamers to your Luminex® posts

HELIX CLASSIC: RSA #1-6, RSB #1-6, RSC #1-5  
STANDARD CLASSIC: RUA #1-6, RUB #1-6.



Dentatus AB  
www.dentatus.com  
SE: +46 8 546 509 00  
info@dentatus.se  
US: +1 212 481 1010  
dentatus@dentatus.com

