

# LUSCENT-**R** ANCHORS™

TRANSLUCENT ENDODONTIC POSTS

## Next Generation Fiber Posts with Increased Radiopacity

### The Advantages

- Radiopaque
- Superior Esthetics
- Dentin-like Elasticity
- Biocompatible
- Priced Right
- Expedited Procedure
- Atraumatic Retrieval
- Swedish Quality

Without compromise, **LUSCENT-R** offers all the benefits expected from brand leading, light-transmitting fiber posts, now with higher radiopacity levels for clear & crisp x-rays.

**Luscent – Twin Luscent – LUSCENT-R**



# LUSCENT-R™ Anchors

Luscent Anchors are a family of fiber-glass, clear resin posts that are designed to refract and transmit natural tooth colors for esthetic post-and-core foundations. **LUSCENT-R** anchors are radiopaque and clearly identifiable on radiographs.

Dentatus Posts have been manufactured in Sweden with micro-precision and unmatched quality for over 70 years.

## EXCELLENT STRENGTH

**LUSCENT-R** is made of high-quality unidirectional fiber rods that are densely distributed and polymerized in a strong resin matrix for superior mechanical strength. The flexural modulus is very similar to that of natural dentin, minimizing risk of root fracture.

## SUPERIOR ADHESION

The adhesion of the epoxy-resin to the silanized fibers in **LUSCENT-R** is a key point for excellent fatigue properties. The texturized surface, in combination with modern adhesive techniques, produces a strong bond to dentin.

## LIGHT TRANSMITTING

The fiber rods effectually transmit light, which allows composite polymerization deep into the canal. When used with dual-cure composites, the additional curing time allows more precise positioning of the post.

## ANATOMICAL DESIGN

The anatomically shaped taper minimizes unnecessary removal of dentin and preserves the natural root structure. **LUSCENT-R** Anchors come in six diameters to provide the optimal choice for the root canal.

## ESTHETIC RESULTS

The neutral shade and translucency of **LUSCENT-R** blend with the tooth, eliminating shadows at the gingival, root and crown interface, as well as, through thin composite restorations.

## ORDERING INFORMATION

### LUSCENT-R STARTER KIT

#### [LUCR-K3]

Contains:

15 posts, 5 each of #4, 5 & 6;

Conical Reamers #4, #5 & #6; and 15 Core-forms

### REFILL PACKS

[LUCR-XXXS10] 10 LUSCENT-R, SIZE 1 & 10 CORE-FORMS


[LUCR-XXS10] 10 LUSCENT-R, SIZE 2 & 10 CORE-FORMS

[LUCR-XS10] 10 LUSCENT-R, SIZE 4 & 10 CORE-FORMS

[LUCR-S10] 10 LUSCENT-R, SIZE 3 & 10 CORE-FORMS

[LUCR-M10] 10 LUSCENT-R, SIZE 5 & 10 CORE-FORMS

[LUCR-L10] 10 LUSCENT-R, SIZE 6 & 10 CORE-FORMS

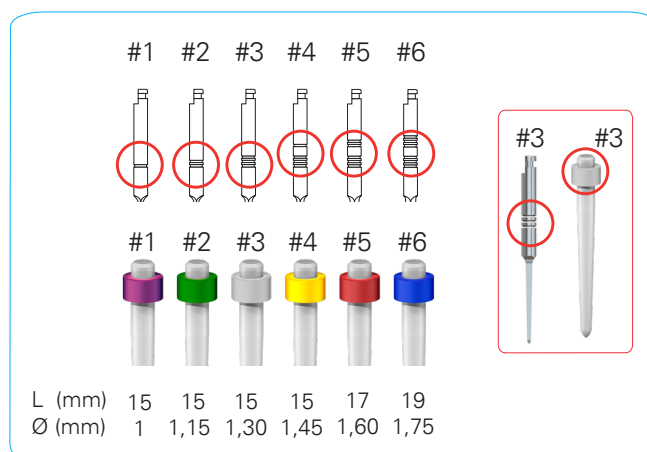

Dentatus AB  
www.dentatus.com  
SE: +46 8 546 509 00  
info@dentatus.se  
US: +1 212 481 1010  
dentatus@dentatus.com

