


DENTATUS COMPLETE IPR SOLUTION

FOR ALIGNER ORTHODONTICS



COMPLETE
IPR
SOLUTION

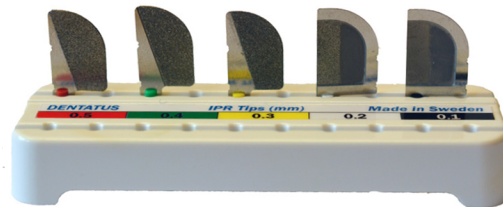

STARTER KIT Item # PDX-P/SK2
includes Profin PDX & Assorted Tip Kit


The Dentatus Complete IPR Solution consists of Interproximal Reduction (IPR) Tips together with the reciprocating Profin PDX contra-angle.

Advantages:

- Significant time savings
- Precise control and accuracy
- Safe to gingival, tongue, lip & cheek
- Comfortable and painless for patient
- Ergonomic design reduces fatigue
- Maintains natural esthetics



Profin PDX with IPR Tip

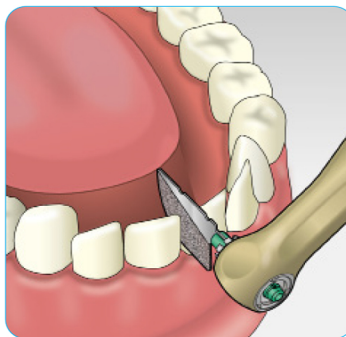

IPR in Align Orthodontics



Clear Aligner Treatment

IPR Tips

Calibrated blade thickness from 0.1 mm to 0.5 mm

Long lasting tips

Autoclavable

CONTROL & PERFORMANCE EXCELLENCE - *Setting a New Standard in IPR*

The Profin PDX handpiece offers unique features such as a reciprocating 1.2mm axial stroke and the ability to “lock” tips, ensuring they are firmly secured in 1 of 8 positions.

SAFETY & CONTROL

With Profin IPR, calibrated tips ensure treatment plan compliance & diminish the risk of excessive stripping. Profin allows for concise treatment with quick and efficient modification.

PERFORMANCE MADE EASIER

Providers can achieve results in 1/6th of the time as compared to manual hand instruments. Profin IPR results in more patient comfort during IPR, higher case acceptance and treatment success.

ESTHETICS

Profin IPR avoids rounding and aids in re-establishing natural embrasures of the tooth. This approach also minimizes scratching & staining and provides a more esthetic result.

STARTER KIT

PDX-P/SK2 includes Profin PDX & Assorted Tip Kit



**Profin
PDX**



**Assorted Tip Kit
LTA-P/SK2**

GETTING STARTED

The easiest way to get started is to purchase the PDX IPR Starter Kit (**PDX-P/SK2**), Included in this kit is the PDX contra-angle, 5 double sided IPR tips and a stand.

If you already own a Profin, consider an Assorted IPR Tip Kit (**LTA- P/SK2**) containing 5 double sided tips 0.1 - 0.5mm or choose from the Refill IPR Tips 2/pk.

REFILL PACKS

IPR tips are color-coded strictly related to blade thickness. This is to ensure time and performance excellence. Grit sizes are noted below.



DOUBLE SIDED IPR Tips

LTA-P1/2: Black

Blade thickness 0.1mm. Grit size 15µ.

LTA-P2/2: White

Blade thickness 0.2mm. Grit size 40µ.

LTA-P3/2: Yellow

Blade thickness 0.3mm. Grit size 50µ.

LTA-P4/2: Green

Blade thickness 0.4mm. Grit size 50µ.

LTA-P5/2: Red

Blade thickness 0.5mm. Grit size 50µ.

LTA-P3/2

SINGLE SIDED IPR Tips

LTA-P25R/2: Violet – abrasive on the right

Blade thickness 0.25mm. Grit size 50µ.

LTA-P25L/2: Violet – abrasive on the left

Blade thickness 0.25mm. Grit size 50µ.



Dentatus AB
www.dentatus.com
SE: +46 8 546 509 00
info@dentatus.se
US: +1 212 481 1010
dentatus@dentatus.com

MADE IN
SWEDEN